



## In-Check DIAL G16 Technical specifications

**Tradition and Innovation** - Visionary thinking and a fascination with technology have guided us throughout our history to develop innovative products of outstanding reliability. Anticipating trends to improve the quality of life.

# In-Check DIAL G16

## Technical specifications



### In-Check DIAL G16

Dimensions	174 x 58mm - w/o mouthpiece, 213 x 58 - with mouthpiece
Weight	106g
Measurement Range	15 - 120 L/min
Accuracy	±10% or 10 L/min
Repeatability	5 L/min

### Product

In-Check DIAL G16	3109300
-------------------	---------

### Disposable Inspiratory One-Way Valve Mouthpieces


Adult Disposable One-Way Valve Mouthpieces Inspiratory 5x20s	3122061
Adult Disposable One-Way Valve Mouthpieces Inspiratory 200s	3122063
Adult Disposable One-Way Valve Mouthpieces Inspiratory 500s	3122065
Paediatric Adaptors for One-Way Valve Mouthpieces (pack of 20)	3122120

### Adaptors

Various adaptors are available to match inhalers not included on the In-Check DIAL G16. Contact us for further information.

### Resistance

High	Med Low
Med High	Low
Medium	pMDI

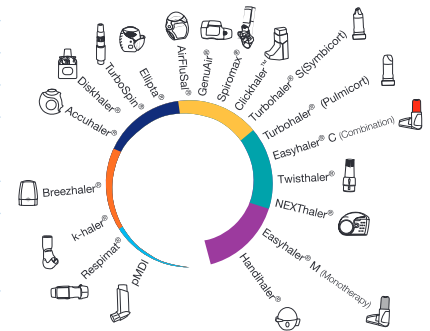

### International alternatives

Handihaler®	AirFluSal®
Easyhaler®	Ellipta®
NEXThaler®	TurboSpin®
Twisthaler®	Diskhaler®
Turbuhaler®	Diskus®
Turbuhaler®, Flexhaler®	Breezhaler®, Aerolizer®, Neohaler®
Clickhaler™	k-haler®
RespiClick®, Spiromax®	Respimat®
Novolizer®, Genuair®,	
Pressair®	

## Product information

The former In-Check DIAL has been revised, now grouping inhalers as a function of their internal airflow resistance, in order to accommodate information on the new devices. There are six resistance groups, five of which relate to dry powder inhalers (DPIs) and one of which to pressurised metered dose inhalers (pMDIs). A guide to identify the resistance group for each inhaler is provided on the In-Check DIAL and on the associated guide. The ability to use In-Check DIAL to assess specific flow rate of a particular inhaler can still be achieved through the use of a specific restrictive adaptor. In-Check DIAL G16 enables the healthcare professional to coach their patients on how to use their inhaler correctly.

- Groups inhalers as a function of their internal airflow resistance, in order to accommodate new inhaler devices
- To coach DPI or pMDI inhaler technique
- Individually calibrated for accuracy
- Easy to clean, by simple washing
- Patients learn to inhale at the correct flow rate
- Inhaler training can be conducted quickly
- Low-cost, simple yet effective method of inhaler technique training

In-Check DIAL G16 is manufactured in accordance with all appropriate EU regulations.

© CLEMENT CLARKE INTERNATIONAL Ltd.

CLEMENT CLARKE INTERNATIONAL Ltd | Edinburgh Way | Harlow | Essex | CM20 2TT, UK  
Phone +44 (0)1279 414969 | Fax +44 (0)1279 456300 | resp@clement-clarke.com | www.clement-clarke.com